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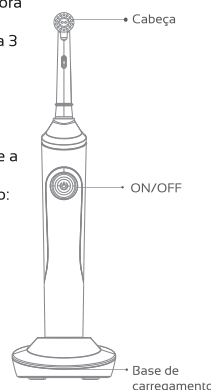
**ATTENTION**

1. Do not put charging base in water or any liquid.
2. Do not recharge, use, or store in fire-hazardous or high-temperature areas.
3. Do not recharge near a strong magnetic field.

**CAUTION**

1. Check if the voltage of the product adapts to your local voltage before using.
2. Replacable items include charging bases, handles, and USB lines. If any item breaks, replace to avoid dangers.
3. Keep out of reach of children. Item is not a toy.
4. Ask doctor before using if you have had an operation in the last two months.
5. Do not use product for other purposes, besides oral care.
6. Clean brush heads with hydrogen peroxide, dicarbonate and bicarbonate.
7. Do not use toothbrush while showering. When not in use, turn off toothbrush to avoid fire.


**USER MANUAL**

1. Align toothbrush head with handle.
2. Push down on head towards handle.
3. Place charging base on a safe and dry surface.
4. Turn on charging base.
5. Place handle on charging base, the indicator light will turn on to show the handle is charging.

**Tip**

Do not dispose of the product when no longer needed. The battery contained inside may be harmful to the environment. Dispose of the product at a recycling center near you.

**SPECS**

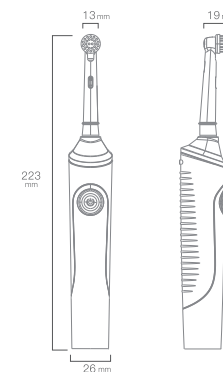- Size 223mm x 27mm
- Colour Azul oceano
- Input power
- Voltage 5V/1A
- Material ABS, TPE, PP

**PROBLEM SOLVING**

SITUATION	REASON	SOLUTION
Vibration reduction	Pushing toothbrush down too hard on tooth. Toothbrush head is secured too close to handle.	Lightly place brush head on tooth. Leave a small gap between brush head and handle.
Short battery life after full charging	Toothbrush is inproperly placed on charging base.	Leave a small gap between brush head and handle. Place handle straight on charging dock and ensure the indicator light is on.
Toothbrush heads break down	Product is used for more than recommended time.	Replace battery.

**BATTERY CHARGE**

1. Charge the toothbrush before using it for the first time.
  2. Charge the toothbrush for approximately 10-12 hours.
  3. It will take approximately 24 hours to fully charge the battery.
  4. Place the toothbrush on the charging base and make sure the indicator light is on.
  5. If the toothbrush has not been recharged for a long time, the indicator light will turn on after 6-10 minutes.
- Normal charge: red indicator  
 Full charge: green indicator  
 Low battery: yellow indicator

**CAUTION**

1. Do not charge in the bathroom. Do not place the charging base near water or shower.
2. Do not break the products or change the battery.
3. Do not change the USB charging line or the charging base.

**BATTERY CHARGE**


**CLEANING THE HEAD**

- Clean the brush head after use and let it air dry.
- Avoid touching the brush to prevent the filaments from unraveling.

**CLEANING THE BODY**

1. Clean the body with a soft cloth.
2. Do not immerse in water.
3. Do not use an acid solvent or water with a temperature higher than 40°.

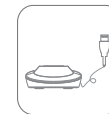**CLEANING THE CHARGING BASE**

1. Remove the base before cleaning to avoid electric shock.
2. Clean dirt and stains with a cloth.

**HEAD REPLACEMENT**

- Replace the head every three months.

**HIGH FREQUENCY**


With a speed of 7600vpm / min, our electric toothbrush will give you excellent and deep cleaning results.

**HEAD**

Toothbrush head: The toothbrush heads are replaceable. The indicator bristles will let you know when it is time to replace the heads (dentists recommend changing the heads every 3 months)

**BUTTONS**


iSwitch: press the switch down, the toothbrush will stop for 0.3 seconds. This prompts you to change the cleaning area. Cleaning time lasts 2 minutes. Charging base: Rechrhes toothbrush.

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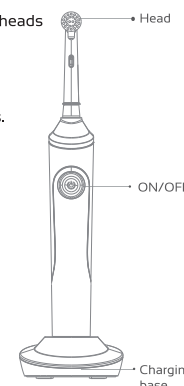## FITTING

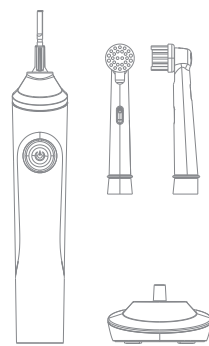Handle	Toothbrush heads	Charging dock
1	2	1

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